

MEDIVATORS® ISA®

Endoscope Reprocessor

Rapid cycle time –
up to 20 minutes

Maximum flexibility
and efficiency

Intuitive handling through
ergonomic design



EVOLUTIONARY PROCESS, REVOLUTIONARY SOLUTION

From more than 20 years, Medivators is dedicated to developing and providing innovative solutions in infection prevention and control worldwide.

As a result of this experience, the **MEDIVATORS® ISA® Endoscope Reprocessor** represents the unique new generation response. It is a rapid, efficient and flexible reprocessing solution designed for endoscopy units demanding the most advanced technology and the highest quality standards.



MEDIVATORS® ISA® ENDOSCOPE REPROCESSOR – THE FAST AND FLEXIBLE REPROCESSING SOLUTION

FAST AND SIMPLE OPERATION FOR IMPROVED EFFICIENCY

- Full cycle in less than 20 minutes with channel monitoring
- Simple, user-friendly handling and operation
- Insertion of endoscopes and operation information through RFID system – helps to avoid errors and saves time




COMPACT AND ERGONOMIC DESIGN FOR MAXIMUM USER COMFORT

- Large and easy to load basin accommodates all types of flexible endoscopes
- Intuitive interface assures easy and fast access to all the cycle information
- Compact size contributes to easy installation, operation and maintenance



ENSURING SAFETY FOR OPERATORS, PATIENTS AND ENDOSCOPES

- Hands free operation and RFID- based technology for complete infection control
- Full traceability through printed and electronic cycle documentation
- Validated process fully compliant to the standards UNI EN ISO 15883-1/4 and UNI CEN ISO/TS 15883/5 (including Austrian and German test soil)

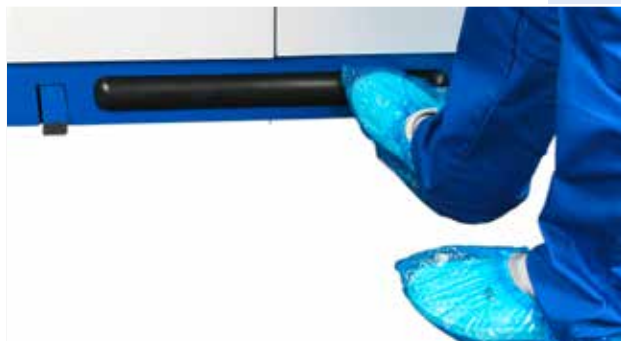

THE RIGHT FORMULA FOR THE BEST PERFORMANCE AND RESULT

- Dedicated room temperature chemicals compatible with all types of flexible endoscopes
- Improved water and chemical consumption
- Designed specifically for use with MEDIVATORS ISA Endoscope Reprocessor



Product Specifications MEDIVATORS® ISA® Endoscope Reprocessor



Specifications

Model	MEDIVATORS ISA Endoscope Reprocessor
Dimensions	700 mm (l) x 1400 mm (h) x 650 mm (d) – monitor included
Reprocessing Time	Up to 20 minutes (full cycle), 12 minutes only (HLD cycle)
Process Capacity	1 endoscope per cycle
High Level Disinfectant/Sterilant	ISASPOR® Single Shot Disinfectant, two parts, Peracetic acid and ISAZONE® Activator. RAPICIDE® PA Single Shot Disinfectant, two part, Peracetic acid.
Detergent	ISACLEAN™ Detergent, multienzymatic detergent INTERCEPT® PLUS Detergent, Non-Enzymatic detergent
Power Requirements	100 V-240 V, 50-60 Hz
Operating Room Temperature	5°-40°C
Weight	75 kg
Humidity	20 -80%

Consumables

Description	Part Number
ISASPOR Single Shot Disinfectant (1x10L Sol. A + B)	500187
ISASPOR Single Shot Disinfectant (1x5L Sol. A + B)	500188
RAPICIDE PA Single Shot Disinfectant (1x5L A + B)	ML02-0158
RAPICIDE PA Single Shot Disinfectant (2x5L AA + BB)	ML02-0115
ISACLEAN™ Detergent (1x10L)	500190
ISACLEAN™ Detergent (2x5L)	500189
INTERCEPT PLUS Detergent (2x5L)	ML02-0145
ISA Sanitary Water Filter Kit	Sanitary: 93250-221 (0,45um + 0,1um)
ISA Colder (CPC) Water Filter Kit	Colder (CPC): 93250-222 (0,45um + 0,1um)
ISA Air Filter	ISA4400
ISA Printer Roll	ISAROTS2

ISA®, ISASPOR®, ISAZONE® and ISACLEAN™ are registered trademarks and trademark of Cantel Medical (Italy) S.r.l.
MEDIVATORS®, RAPICIDE® and INTERCEPT® PLUS are registered trademarks of Medivators Inc.



Cantel Medical (Italy) S.r.l.
Via Laurentina, 169
00071 Pomezia (RM) Italia
Tel.: +39 06 9145399
Fax: +39 06 9146099
www.cantelmedical.it

Medivators Inc.
14605 28th Avenue North
Minneapolis, MN
55447-4822 USA
Tel: +1.763.553.3300
Fax: +1.763.553.3387

Medivators BV
Sourethweg 11
6422 PC Heerlen
The Netherlands
Tel: +31.45.5.471.471
Fax: +31.45.5.429.595

Cantel (UK) Limited
Campfield Road
Shoeburyness, Essex
SS3 9BX,
United Kingdom
Tel: +44 (0) 1702.291878
Fax: +44 (0) 1702.290013

Cantel Medical Asia/Pacific Pte. Ltd.
1A International Business Park
#05-01 Singapore 609933
Tel: +65.6227.9698
Fax: +65.6225.6848

Cantel Medical Devices (China) Co. Ltd.
Unit 804-805, Innov Tower Block A,
Hongmei Road, Xuhui 200233 Shanghai
Tel: +86 21 60161380
Fax: +86 21 61210913

www.cantelmedical.eu

© 2017 Medivators Inc. 50098-882-EN/E