

SURESTORE[®]

Endoscope Storage & Transportation | **Dry & Store**





THE COMPLETE CIRCLE OF PROTECTION

As the global vanguard in infection prevention, only Cantel delivers the Complete Circle of Protection, a full-value, proactive partnership dedicated to helping you remove risk, streamline operational efficiencies and optimise your success.

Cantel is proud to introduce the SURESTORE® Endoscope Storage & Transportation

British designed and manufactured by our experienced engineers and technical staff, the SURESTORE® Endoscope Storage and Transportation solution offers cutting edge long-term 'active' aseptic storage and transportation of endoscopes (including GI and theatre based scopes). Conditioned directly from an endoscope washer disinfector such as RAPIDAER® Endoscope Reprocessor, and completely negating the need for time consuming air drying, the SURESTORE® Endoscope Storage and Transportation solution validates each scope for storage of up to 100 days.

The unique design and technology of the SURESTORE® Endoscope Storage and Transportation solution provides a 'patient ready' packaged endoscope within 10 minutes.

SURESTORE®
Storage & Transportation



RAPIDAER®
Endoscope Reprocessor




CLEANSCOPE®
Transport & Short
Term Storage System



Conditioned
Packaged
Endoscope



Hospital Transportation

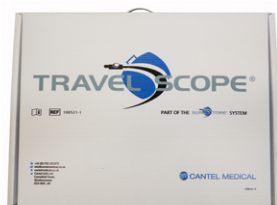

National Transportation



Regional Transportation



DRY & STORE

Bacteria pose significant risk to endoscopes during endoscope storage. Cantel's transport, drying and storage solutions are designed to protect valuable inventory, reduce cross contamination touchpoints, eliminate moisture in the endoscope channels and control humidity. Higher humidity or moisture in endoscope channels is known to aid bacterial growth.

SURESTORE®

Endoscope Storage & Transportation | **Dry & Store**

Convenient touch screen operation

Visible cycle progress information and independent monitoring

Contoured tray for safe & secure endoscope storage (compatible with CLEANASCOPE and the TRANSCOPE/STORASCOPE lockable carts)

Compatible with the CLEANASCOPE sterile packed liner system to reduce the potential for cross-contamination

Lockable drawer with keyboard to input new information & hidden compact PC for space saving

Fully traceable print cycle data

Large lockable storage drawer & hidden safe chemical storage

Independently lockable castors



SURESTORE® Endoscope Storage and Transportation solution enables scopes to be easily prepared in advance for emergency situations, out of hours, weekend, and off-site endoscopy work. The scopes can be transported to various departments within a hospital, or to secondary sites, using the TRAVELSCOPE® Transportation Case, or the TRANSCOPE® & STORASCOPE® Lockable Carts.

- Validated storage for 100 days
- Independent monitoring system
- Reduced scope wear and tear
- Flexible transportation solution
- Fully traceable cycle data
- Quick and simple to use
- Secure and robust storage
- Space saving storage solution
- Independently validated by Biotech Germande
- Tested with multiple endoscope brands

Product Specifications SURESTORE® Endoscope Storage & Transportation



Specification

Model	SURESTORE® Endoscope Storage & Transportation
Conditioning time (start to finish)	Approximately 8 minutes (lumen scopes) Approximately 2 minutes (non-lumen scopes)
Dimensions	530mm (w) x 740mm (d) x 1360mm (h)
Power requirements	13 Amp 220 - 240V
Operating requirements	Completely mobile on 4 independently lockable castors
Installation space	Unit does not require any additional ventilation space and can therefore be positioned against a wall or between other equipment
Validation	Biotech Germande and NFS 98-030

SURESTORE® Endoscope Storage & Transportation ordering information

Part number	Description	Units per box
100600	1 X SURESTORE® Storage/Transportation Unit 12 X SURESTORE® Supportive Trays 1 X 5 Litre N-SURE® Solution 50 X SURESTORE® Outer Pouches 45 X SURESTORE® Sterile Collection Pads 50 X CLEANASCOPE® Liners 1 X Roll of 300 Tracking Labels	1
100602	N-SURE® Solution 5L	1
100710	N-SURE® Zero Packs (non-lumen scopes)	50
100604	SURESTORE® Outer Pouches	50
100728	SURESTORE® Sterile Collection Pads	225
100606	SURESTORE® Supportive Tray	1
700344	SURESTORE® Tracking Labels	300
200018	SURESTORE® HEPA Air Filter	1
103580	CLEANASCOPE® Sterile Liner Set (red and green covers with base liner)	200
100521	TRAVELSCOPE® Transportation Case	1
103551	TRANSCOPE® Lockable Cart 6 Tray (contoured trays)	1
103543	TRANSCOPE® Lockable Cart 6 Tray (standard trays)	1
103544	TRANSCOPE® Lockable Cart 6 Tray (standard trays with lids)	1
103553	STORASCOPE® Lockable Cart 10 Tray (contoured trays)	1
103547	STORASCOPE® Lockable Cart 10 Tray (standard trays)	1
103548	STORASCOPE® Lockable Cart 10 Tray (standard trays with lids)	1

SURESTORE®, CLEANASCOPE®, TRAVELSCOPE®, TRANSCOPE®, STORASCOPE®, and N-SURE® are registered trademarks of Cantel (UK) Limited.



Cantel (UK) Limited
Shoeburyness
Campfield Road
Essex, SS3 9BX
United Kingdom
Tel: +44 (0) 1702 291878
Fax: +44 (0) 1702 290013

Medivators BV
Sourethweg 11
6422 PC Heerlen
The Netherlands
Tel: +31.45.5.471.471
Fax: +31.45.5.429.595

Medivators Inc.
14605 28th Avenue North
Minneapolis, MN
55447-4822 USA
Tel: +1.763.553.3300
Fax: +1.763.553.3387

**TO PLACE
AN ORDER**

Call: 01702 291878
Fax: 01702 290013
Email: orders@cantelmedical.co.uk